


Priority Principle™: Sciatica




Priority	Principle	Frequency	Time
0	Pain (PRN)	(PRN)	
2nd	Inflammation	50, 5-250 Hz	DOSE
3rd	Spasms	1000 Hz PRN (treat bilaterally if applicable)	Pontinen's Principle
4th	(a) Pain (Systemic)	NRT (Dermatome) 1000, or 1000-3000, 3000, 5000 Hz	4 minutes scanning along the affected dermatome
	OR		
(b)	Pain (Local)	1000, 3000 or 5000 Hz	3-5 minutes

# of treatments					

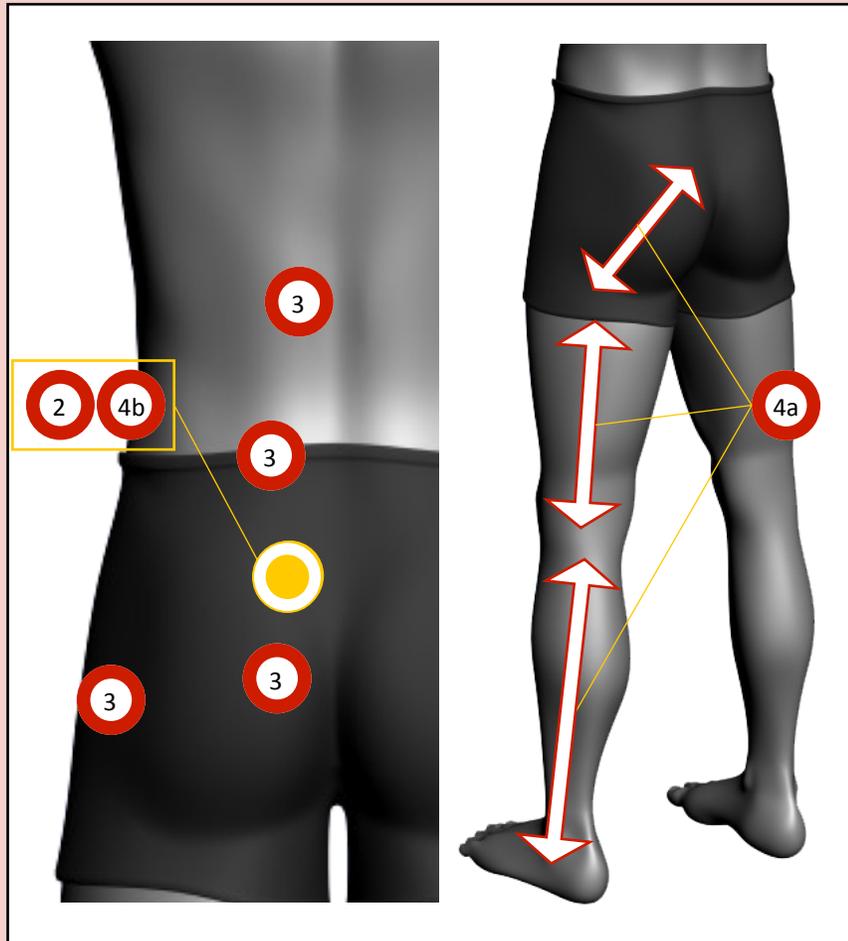



Priority Principle™: Sciatica

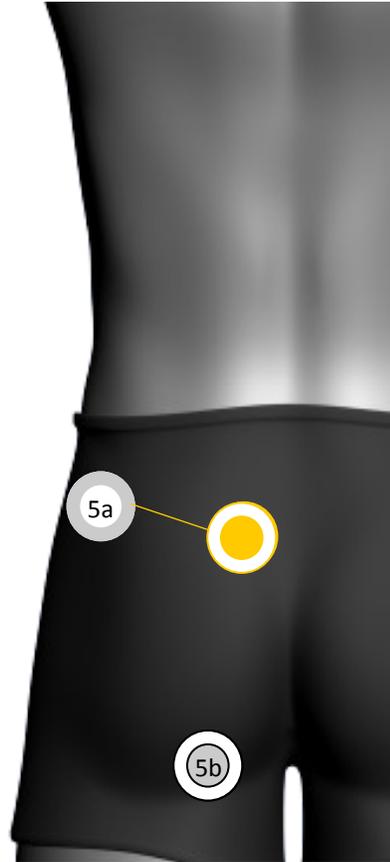

Priority	Principle	Frequency	Time
5th	(a) Tissue Repair (Primary)	250, 5-250 or 5-500 Hz	DOSE
	and		
	(b) Tissue Repair (Secondary)	50 Hz PHT (Femoral Artery)	5 minutes

# of treatments					

Priority Principle™: Sciatica

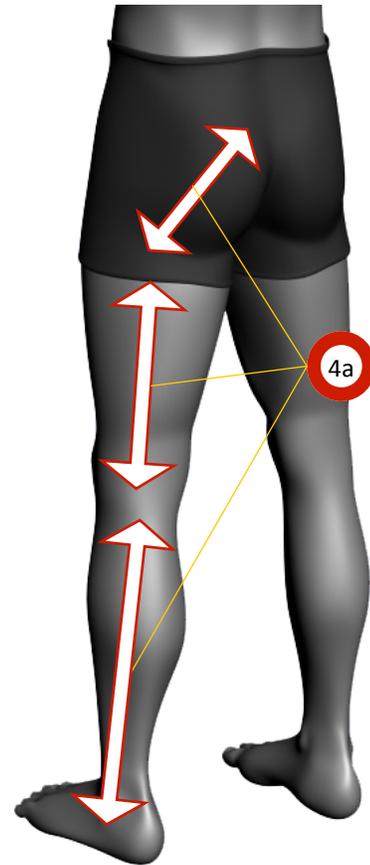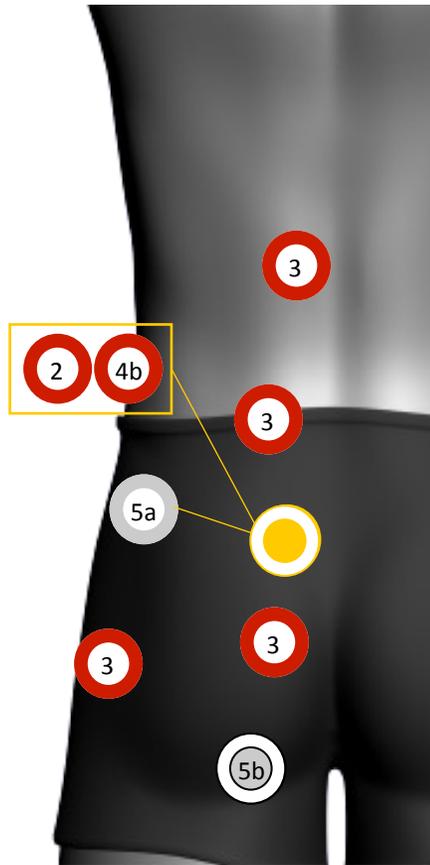

Priority	Principle	Frequency	Time
0	Pain (PRN)	(PRN)	
2nd	Inflammation	50, 5-250 Hz	DOSE
3rd	Spasms	1000 Hz PRN (treat bilaterally if applicable)	Pontinen's Principle
4th (a)	Pain (Systemic)	Lumbar Sympathetic Plexus Laser Block	See Method (to be done bilaterally)
OR			
(b)	Pain (Local)	1000, 3000 or 5000 Hz	3-5 minutes
5th (a)	Tissue Repair (Primary)	500-1000 Hz	DOSE
and			
(b)	Tissue Repair (Secondary)	50 Hz PHT (Femoral Artery)	5 minutes

